



Hyagen



Stick

Welcome to the world of Hyagen® Stick

Since the 1950s, the first attempts have been made with electric current in cosmetic treatment. From the 90s onwards, the methods were more and more refined.

The electro-physical cosmetic treatment has been further developed into an effective treatment method in the anti-aging field. Leading scientists believe that electro-physical cosmetics can bridge the gap between traditional and traditional cosmetics in terms of mode of action and results. cosmetics and cosmetic surgery.





What is electrocosmetics?

The human body is a "liquid conductor of electricity". All body fluids contain ions, which are also permanently formed (e.g. by dissolving salts). If a small electrical voltage is applied to the skin via an anode plate, the electrically charged particles (ions) already present in the body begin to migrate to the anode plate, depending on their electric charge.

Directly below the surface of the skin, movement is created by the applied tension.

Care products are optimally incorporated into the skin through this permanent movement of the ions. The effectiveness of an order from Care products with the Hyagen® Stick exceed those of a hand application many times over.

Furthermore, the ion movements calm nerve receptors and improve blood circulation and lymph flow. Hyagen® is the care product perfectly matched to the Hyagen® Stick.



Function modes:

Active ingredient absorption / massage

Clinical studies have shown that a maximum of 15% of the active ingredients present in creams can penetrate the skin barrier when applied by hand. If the Hyagen® Gel, specially developed for the Hyagen Stick, is now used and worked into the skin with the Hyagen Stick by means of activated "Active Ingredient Absorption / Balance" mode, up to 80% of the active ingredients can pass through the skin barrier in order to be able to fully develop the effect.

Hyagen is a gel specially developed for the Hyagen Stick, which has excellent care properties and is enriched with special active ingredients developed for electro-physical cosmetic treatment.

The result is cells that absorb more moisture and thus restore freshness and elasticity to the skin while at the same time reducing wrinkles.




Application

By pressing the switch "Up arrow" turn on the device. The treatment head for the eye area is activated.

When the "arrow down" switch is pressed, the treatment head for the face is activated.

To switch off the device, move the switch to the middle position (neutral position)

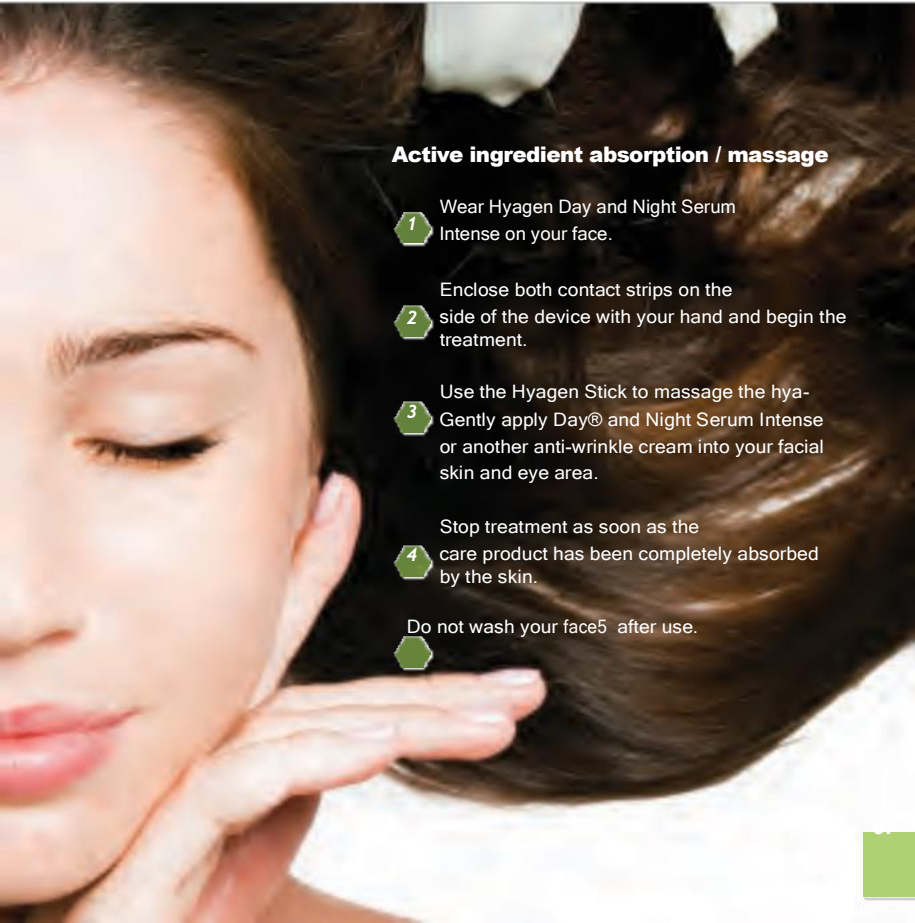

Active ingredient absorption / massage

1

Wear Hyagen Day and Night Serum Intense on your face.

2

Enclose both contact strips on the side of the device with your hand and begin the treatment.

3


Use the Hyagen Stick to massage the hya- Gently apply Day® and Night Serum Intense or another anti-wrinkle cream into your facial skin and eye area.

4

Stop treatment as soon as the care product has been completely absorbed by the skin.

5

Do not wash your face⁵ after use.

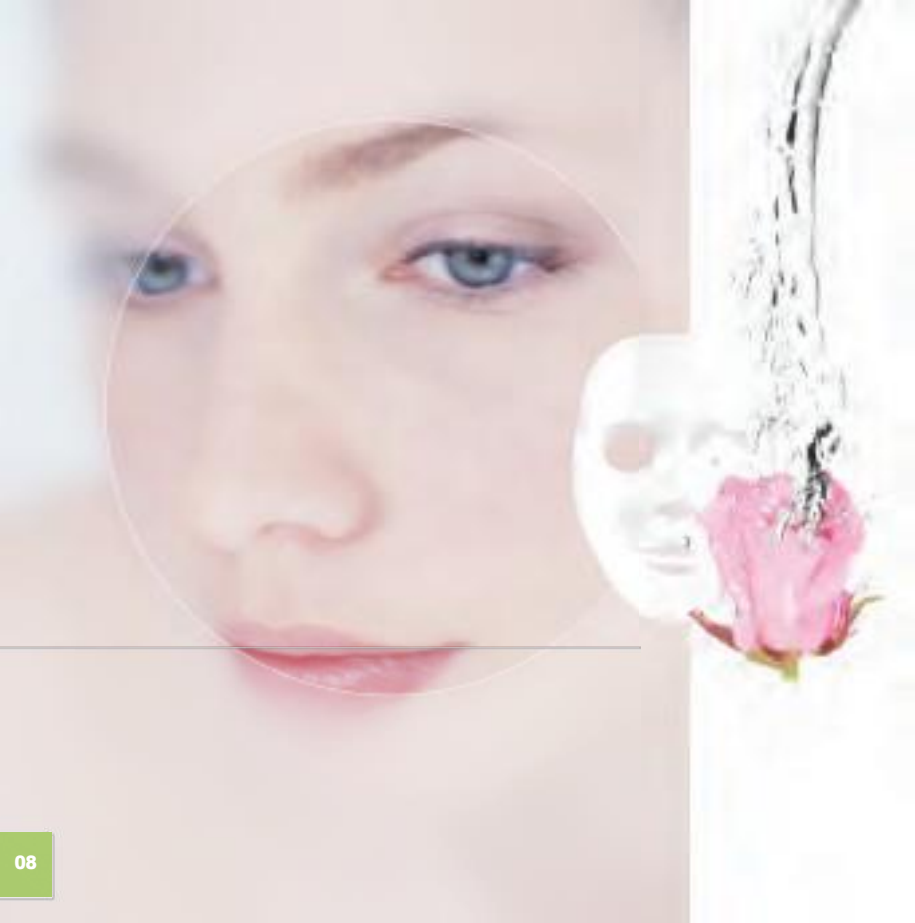

Skin

Neutral skin 1 time a day; 5 – 10 minutes

Sensitive / Oily Skin 1 time every 3-4 days; 5 – 10

minutes

Extremely sensitive /
Extremely oily skin 1 time a week, maximum 5 minutes

Aged skin 1 time a day; 10 - 15 minutes

Eyes The skin around the eyes is very sensitive.
Remain on a point for a maximum of 2
seconds.

Hints

Always touch the two lateral ion contact strips of the device when using it.




Hyagen
Parfums de haut de gamme
Intense



100% natural

Hyagen® Day and Night Serum Intense

Hyagen® is a gel specially developed for the Hyagen® Stick, which has excellent care properties and is enriched with special active ingredients developed for electro-physical cosmetic treatment.

Hyaluronic acid is naturally stored in the skin as well as in other areas of our body. The substance is produced by the human body itself and is characterized by up to 1000 times the storage of water in relation to its own weight, thus ensuring optimal hydration and firmness of the skin.

With increasing age, the production of hyaluronic acid in the body decreases continuously. The skin loses density, elasticity and suppleness. wrinkles arise. External influences, such as strong sunlight or pollution from the environment and dust, can promote the appearance of skin fatigue and accelerate the aging process.

Due to the increased moisture absorption, Hyagen® has a particularly skin-tightening effect. Hyagen® - Day and Night Serum is 100% natural and contains no animal ingredients. Hyaluronic acid is found in all cells of the human body; allergic reactions are therefore extremely rare.

Thanks to biotechnological manufacturing processes, hyaluronic acid of animal origin is dispensed with.

Maintenance

- Regular maintenance and cleaning ensures optimal results and extends the service life of the device.
- Store Hyagen Stick in a safe, cool and dry place when the device is not in use.
- In order not to compromise the safety, reliability and performance of the device, no extensions, readjustments, modifications or repairs may be carried out.
- Extensions, new hires, changes or repairs may only be carried out by persons authorised by Canbolat Vertriebs GmbH.
- Before cleaning the appliance, make sure that it is switched off and unplugged.
- For cleaning, use a slightly damp and soft cloth.
- Do not use aggressive substances such as scouring milk or other washing-up agents, as these can attack the surface.
- Do not repair the device. It does not contain any parts that can be serviced by the user.
- Remove cosmetics or other residues from the Hyagen Stick after each use.
- Clean the appliance regularly with a clean cloth.
- Do not use water to clean the device

Warnings

- This device is intended for domestic use.
- The device is not waterproof - do not contact it with water.
- Store the device in a dry environment.
- Make sure that the device is not exposed to strong vibrations
- Do not open the device or attempt to repair it.
- People with metallic implants, pacemakers or pregnant women should consult a doctor before use.
- The device must not be used on wounds, eczema or tumors. Keep the device out of the reach of children.
- Remove piercings, earrings or other metallic objects before use.
- Do not pass the switched on device over the eyes.
Do not use the device if you have allergies to metal.
- Never place or use the appliance in a wet or damp environment.
- Do not use the device under water, e.g. in the shower.
- Never bring the appliance into direct contact with fire, gas, oxygen or hot objects such as stovetops.

- Handle the device carefully so that it cannot fall or otherwise be damaged.
- If there are any problems with the device, take it in for repair immediately. Don't try to fix it yourself.
- Unplug the AC adapter when not in use.
- This device is not intended for use by persons (including children) with reduced physical, sensory or mental abilities or with a lack of experience and/or knowledge, unless they are supervised by a person responsible for your safety and have been instructed on how to use the device.
- Do not carry the device by the cable or use the cable as a handle.
When you exit the application, set all switches to Off/Off/0 and unplug the power supply.
- Only use extension cords if you have previously made sure that they are in good condition.
If the power cord of the device is damaged, it must only be replaced by an authorized service center. Never let children play with the
- packaging material. There is a risk of suffocation.
- The Hyagen Stick may only be used for its original intended purpose.

Disposal

Power tools, accessories and packaging are to be recycled in an environmentally friendly way.



Do not throw power tools in the household waste! According to the European Directive 2002/96/EEC on waste electrical and electronic equipment and its transposition into national law, power tools that are no longer usable must be

are collected separately and recycled in an environmentally sound manner.

Within the EU, this symbol indicates that this product must not be disposed of in household waste. WEEE contains valuable recyclable materials that should be recycled and so as not to harm the environment or human health through uncontrolled waste disposal.

Therefore, please dispose of old appliances via suitable collection systems or send the appliance to the place where you bought it for disposal. This will then send the device for material recycling.

Accumulators and Batteries

Do not throw rechargeable batteries in household waste, water or fire. Batteries should be collected, recycled or disposed of in an environmentally friendly way.

Only for EU countries

In accordance with Directive 91/157/EEC, defect or used rechargeable batteries/batteries can be handed in directly to Canbolat Vertriebsgesellschaft mbH, Gneisenastr. 10-11, 97074 Würzburg.

Specification of spare parts to be used

The device is delivered without spare parts. If necessary, please contact our service at +49 931/90803000.

Technical data

Product: Hyagen Stick

Dimensions: 15.5 x 4 x 3.5
cm

Rotation Per Minute: 12000 Power

Consumption: 0.2W

